
medisana®

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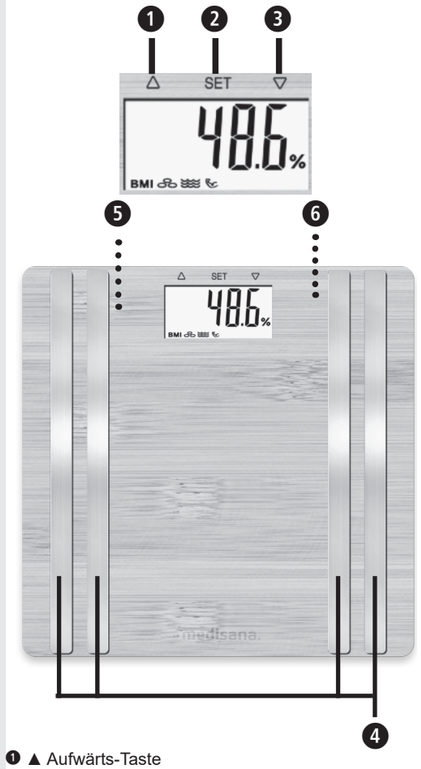
medisana®

GB Instruction manual
Body analysis scale BS 550 / BS 552

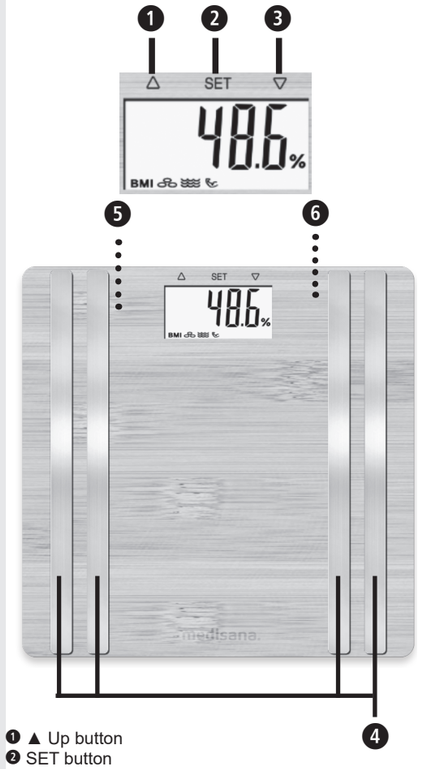❶ ▲ Up button
❷ SET button
❸ ▼ Down button
❹ Electrodes
❺ Battery compartment (on the underside)
❻ Switch for weight units (on the underside) for kilograms (kg), pounds (lb), stones (st)

Legend

IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes give you useful additional information on the installation or operation.

LOT number

Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB Important Instructions Retain for future use

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- This scale is not suitable for people with pacemakers or other medical implants. The indication of body fat may be inaccurate in the event of diabetes or other medical/physical conditions. The same applies to highly trained sportsmen.
- This device is not suitable for pregnant women!
- Any treatment or diet for overweight or underweight requires qualified advice from appropriate specialists (doctor, dietician). The values determined with the scale can be of help.
- Risk of tilting! Never place the scales on an uneven surface. Never step onto just one corner of the scale. While weighing, never stand on one side or on the edge of the scale.
- Danger of slipping! Never step onto the scale with wet feet. Never step onto the scale while wearing socks.
- Use the device only for its intended purpose, as directed by the user instructions. In the event of misuse, all guarantees become null and void.
- The scale is manufactured for home use. It is not intended for commercial use in hospitals or other medical facilities.
- Do not use the device if it does not work properly, if it has been dropped, dropped into water or damaged.
- The scale has a measuring range up to 180 kg / 396 lb / 28 st 4 lb. Do not overload the scale.
- Place the scale on a firm level surface. Soft, uneven surfaces are unsuitable for measurement and lead to incorrect results.
- Place the scale in a location where neither extreme temperatures nor extreme humidity are present.
- Keep the scale away from water.
- Please handle the scale with care. Avoid shocks and vibrations of the scale. Do not drop it.
- Stand carefully on the scale. Do not hop or jump on the tread surface, as this may damage the weighing mechanism.
- Do not attempt to disassemble the device, otherwise the warranty will become void. This device does not contain any parts that can be serviced or replaced by the user.
- In the event of malfunctions, do not proceed with any unauthorised repair of the device. Otherwise, this will render any warranty claims void. Have repairs performed only by an authorised service centre.

Battery safety instructions
- Do not take the batteries apart!
- If required, clean the battery and device contacts before insertion!
- Remove spent batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! In the event of contact with battery acid, immediately rinse the affected area with plenty of clean water and consult a doctor immediately!
- If a battery is swallowed, consult a doctor immediately!
- Insert the batteries correctly, paying attention to the polarity!
- Keep the battery compartment closed properly!
- If the device is not used for an extended period of time, remove the batteries!
- Keep the batteries out of the reach of children!
- Do not re-charge the batteries! There is a risk of explosion!
- Do not short-circuit the device! There is a risk of explosion!
- Do not throw into a fire! There is a risk of explosion!
- Store unused batteries in the packaging and not near metal objects to avoid a short circuit!

Scope of delivery
Please check first whether the unit is complete and does not show any damage. If in doubt, do not operate the unit and contact your dealer or service centre.
The scope of delivery includes:
- 1 **medisana** Body analysis scale **BS 550** / **BS 552**
- 3 Batteries (Type AAA, 1.5 V)
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Requirements for correct measurement results
- Always weigh yourself in the same conditions.
- Measure at the same time of day.
- Stand upright and calm.
- It is best to take the measurement in the morning, after showering or bathing, with dry feet.

With the results of the measurement, when carried out consistently, a reliable control of the development of the body weight is possible.

Initialise scale
Make sure that the scale is on firm and level ground. If you want to weigh yourself, but have moved the scale immediately before, you must first initialise the scale. Briefly press the middle of the scale surface with your foot. The display will show “0.0”. When the scale has turned off, it is ready for the “step-on” function. If the scale has not been moved beforehand, this process is omitted.

Connect with the VitaDock+ App
Download the VitaDock+ App free from the AppStore or the Google Play Store. Activate Bluetooth on your iOS or Android device. Start the App. Register now with your VitaDock online account. If you don’t have your own account yet, you can register free now at VitaDock Online. Tap the Bluetooth icon to go to the device selection. Scroll down, until the body scale is displayed in the device list. Tap the arrow next to
the illustration to start the setup process.

Simple weighing
- 1. Stand on the scale and stand still. The device switches on automatically and the display briefly shows “0.0”.
- 2. Your weight will be measured and displayed.
- 3. Step down off the scale. The scale switches off automatically after approx. 15 seconds.

Body analysis
The scale can store the data of 8 persons: gender, age, height.
- 1. Briefly press the middle of the scale surface with your foot. The display will show “0.0”.
- 2. Long press the SET button ❷. The memory location flashes in the display.
- 3. Now you can select the desired memory location for your user profile by pressing button ❶ or button ❸.
- 4. Press the SET button ❷ to save the setting. The symbol for the gender starts flashing in the display.
- 5. Set your gender by pressing button ❶ or button ❸ to select the corresponding symbol in the display.
- 6. Press the SET button ❷ to save the setting. The preset height begins to flash in the display.
- 7. Now set your height with the button ❶ or key ❸.
- 8. Press the SET button ❷ to save the setting. The preset age begins to flash in the display.
- 9. Use the button ❶ or the button ❸ to set your age.
- 10. Press the SET button ❷ to save the setting.
- 11. The configuration process is now complete. The display shows “0.0”. Stand on the scales barefoot and stand still. Place your feet on the electrodes ❹. Firstly, the weight is displayed. Thereafter, the measured values for body fat, water content, muscle percentage, bone weight, BMI and calorie requirement (BMR) are displayed one after the other.

Change the weight unit
Briefly press the middle of the scale surface. The display shows «0.0». You can now select the desired units (kg - lb - st) with the button on the back of the scale.

Delete user data
Briefly press the middle of the scale surface with your foot. The display will show “0.0”. Press the down button ❸ for 3 seconds until the display shows ‘dEL’. Confirm deleting the user data by pressing the SET button ❷. The display ‘-----’ flashes twice: All saved data have been deleted.

Error messages
O-Ld = The scale is overloaded.
LO = The battery is weak and needs to be replaced.

Cleaning and care
Never use aggressive cleaning agents or stiff brushes. Clean the scale only with a soft, slightly damp cloth. Never use harsh detergents or alcohol. No water should enter the device.
Only use the device again when it is completely dry.

Technical data
Name and model: **medisana** Body analysis scale **BS 550** / **BS 552**
Power supply: 4.5 V (3x 1.5 V, AAA, Micro LR03)
Measuring range: up to 180 kg/ 396 lb / 28 st 4 lb
Division: 0.1 kg / 0.2 lb / 0.2 lb
Dimensions: 300 x 300 x 35 mm
Weight: 1.35 kg
Item No.: 40550 / 40552
EAN No.: 40 15588 40550 1 / 40 15588 40552 5
Frequency band: 2400 ~ 2478 MHz
Maximum radio-frequency power transmitted in the frequency band: 8 dBm

Guidelines and standards
We hereby declare that this device is in compliance with the essential requirements of the European RED 2014/53/EU directive. The complete EU declaration of conformity can be requested from medisana GmbH, Carl-Schurz-Str. 2, 41460 Neuss, Germany.

Disposal
This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms
Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.
The following warranty terms apply:

1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH,
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under https://docs.medisana.com/40552.
Information about service can be found here: www.medisana.com/servicepartners

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